

FilmArray[®] 2.0 System

The Fastest Way to Better Results

Comprehensive
Test Panels



Intuitive
Step-by-Step
Software



FilmArray
Instrument

1
Hour

Results in about 1 Hour

FilmArray: Fastest Way to Better Results

Bench Top System – This all-in-one system integrates sample preparation, amplification, detection and analysis into one instrument that processes sample to results in about one hour.


Efficient – With just two minutes of hands on time, this closed system does all of the work so you can start the run and return for results.

Comprehensive Test Panels – A system that can identify dozens of the most lethal viruses, bacteria and parasite including emerging infectious diseases. This results in a revolutionary detection system with a lightweight, small-footprint format.

A Small-Footprint for any Mobile Lab

Single Instrument Integration - Reduce the amount of equipment and consumables with the FilmArray System. With its small footprint, it is the ideal instrument for any mobile lab that is on the go.

Freeze-Dried Reagents - All reagents have been freeze-dried and are stable at room temperature.

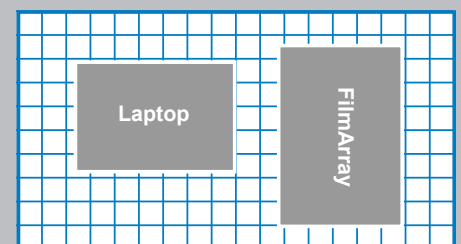

36" Bench Space

Expand with the FilmArray as your Detection Needs Expand

FilmArray has set the standard for pathogen testing. With single database management of up to 8 instruments per computer, FilmArray can test up to 176 samples a day. Enabling large scale test ordering, quicker turnaround times and increased accuracy by minimizing manual data entry.

- High throughput
- Scalable
- Random access performance
- LIS capable
- Single database management



8
FilmArrays

Network up to 8 FilmArray Instruments per computer to test up to 176 samples per day

Intuitive Step-by-Step Software

The FilmArray Software includes a step-by-step wizard that walks the operator through loading the pouch into the instrument and demonstrates how to perform a run. When the run is complete, the advanced software automatically determines and reports results in a simple and easy to read format.

Note: The FilmArray 2.0 software has been designed specifically for the FilmArray 2.0 instrument and laptop and is not backwards compatible with older FilmArray Systems.



A Mobile Lab in Hard Case

If you are testing samples near the hot zone, in the field or out of a mobile command center, then our hard case is what you need. This optional case provides the durability and ruggedness to keep a FilmArray System protected and ready for use. All of the components fit into this custom built hard case.



FilmArray 2.0 Instrument and Components

Part No.	Description
FLM2-ASY-0001	FilmArray 2.0 Instrument Includes: US and European Power Cords, Ethernet Cable, Pouch Loading Station, FilmArray 2.0 Software Installation and Information CDs.
COMP-SUB-0030	FilmArray 2.0 Laptop Computer with US Power Cord Includes: Microsoft Windows 7 Professional 64-bit Edition with FilmArray System Software, Intel Core i5, 4 GB RAM Memory and a 500 GB HDD.
FLM2-ASY-0101	FilmArray 2.0 Laptop Peripherals Kit Includes: Laptop Manual Addendum to the System Manual, Mouse for computer and a programmed Barcode Reader.
FLM2-ASY-0100	FilmArray 2.0 Network Kit Includes: 16-port, 10/100 Mbps Network Switch, Ethernet Cable and USB Network Adapter
NGDS-CAS-0001	FilmArray 2.0 Rugged Hard Sided Case ** optional case, not included with standard FilmArray package